

# Navigator™ HD


Ureteral Access Sheath

Boston  
Scientific

## TECHNICAL DATA SHEET



### Dilator Tip Flexibility

The Navigator HD Ureteral Access Sheath retains the characteristics of the flexible tip of the current Navigator Ureteral Access Sheath, while still providing a stiffer shaft for ease of placement.

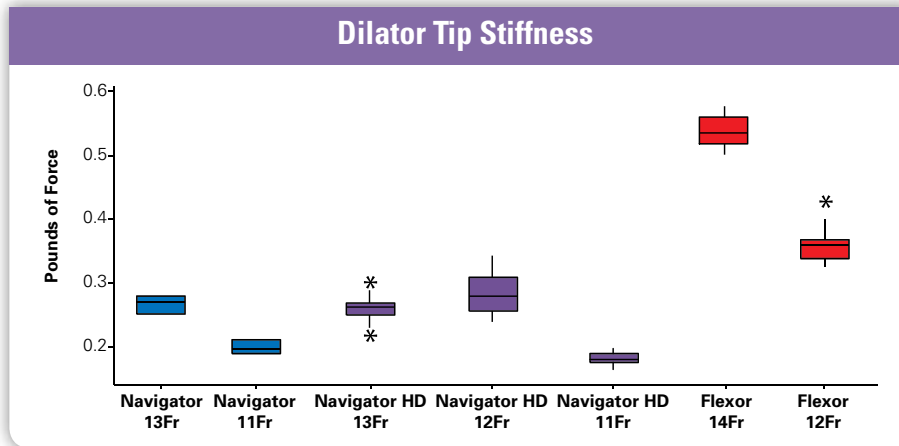

### Improved Coating


The hydrophilic coating on both Navigator HD Ureteral Access Sheath and dilator show an improved coefficient of friction and an improved durability over Navigator Ureteral Access Sheath and Cook Flexor™ Ureteral Access Sheath after undergoing 20 cycles of placement through a tight simulated stricture.

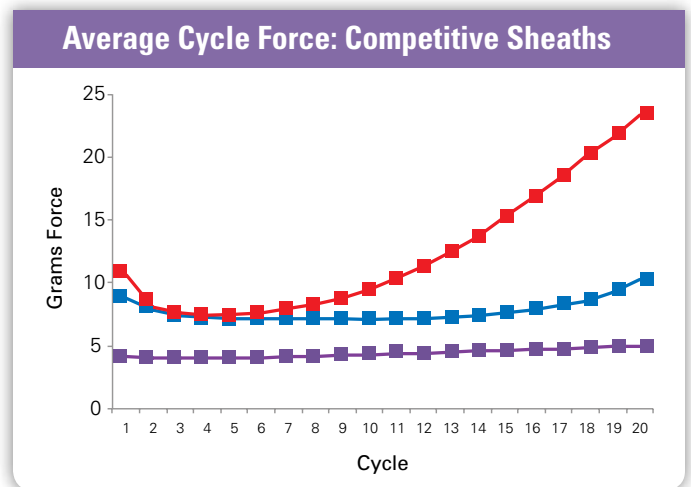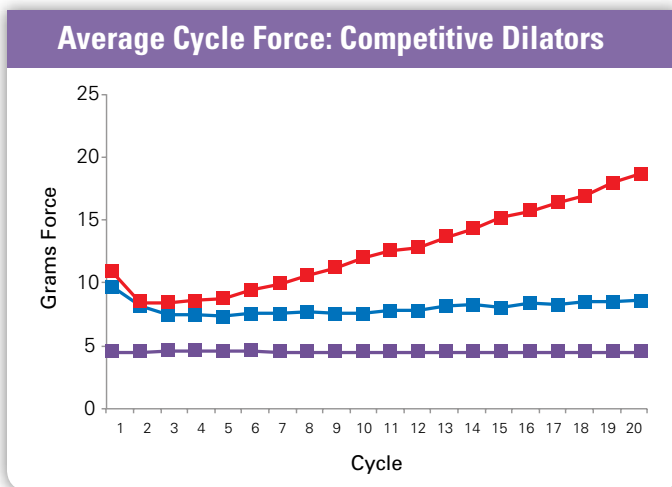

**Color Key:** ■ Navigator Sheath ■ Navigator HD Sheath ■ Flexor Sheath

## Lower Insertion Force

The Navigator™ HD Ureteral Access Sheath requires less force to insert through a simulated stricture when compared to Navigator Ureteral Access Sheath and to Cook Flexor™ Ureteral Access Sheath, indicating a smoother transition from guidewire to dilator and from dilator to sheath.

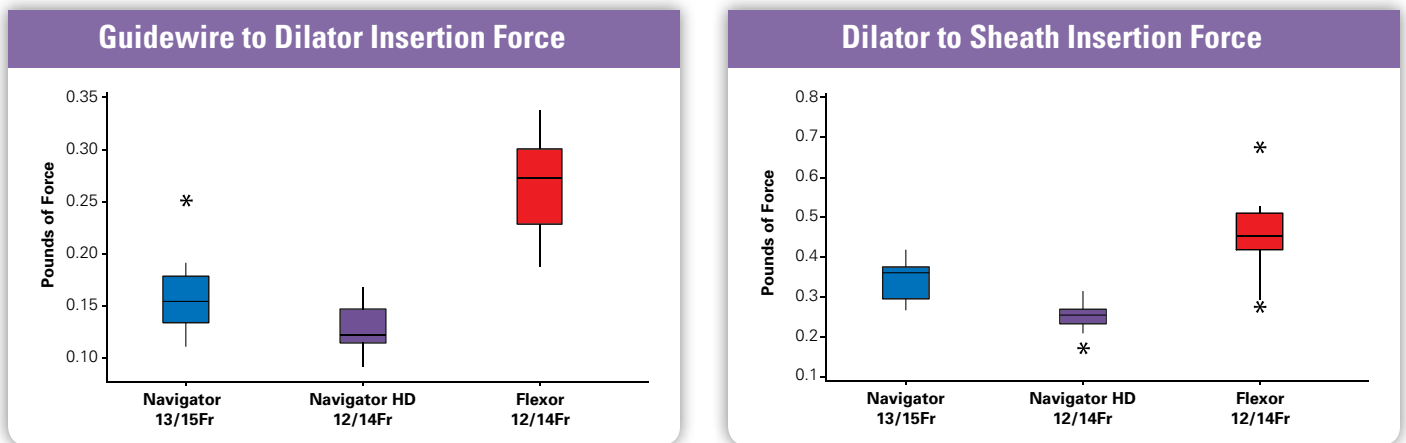


**Color Key:** ■ Navigator Sheath ■ Navigator HD Sheath ■ Flexor Sheath

## Hub Design

### One handed latch engagement:

The Navigator HD Sheath can be successfully assembled with the dilator in almost any rotation relative to the sheath hub vs Cook Flexor Sheath can only be successfully assembled when the dilator is aligned rotationally with the sheath hub.

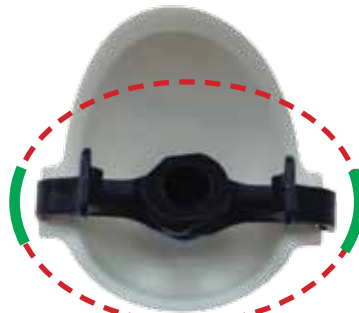

Cook Flexor Sheath



Navigator HD Sheath

**Boston Scientific**  
Advancing science for life™

Boston Scientific Corporation  
One Boston Scientific Place  
Natick, MA 01760-1537  
www.bostonscientificstone.com

**Ordering Information**  
**1.888.272.1001**

©2013 Boston Scientific Corporation or its affiliates. All rights reserved.

URO-132307-AA XM 1/13

Data on file at Boston Scientific Corporation and will be made available upon request.

Navigator is a registered trademark of Boston Scientific Corporation or its affiliates. All other trademarks are property of their respective owners.

**Caution:** Federal Law (USA) restricts this device to sale by or on the order of a physician. Refer to package insert provided with the product for complete Indications for Use, Contraindications, Warnings, Precautions, Adverse Events, and Instructions prior to using this product.